

MStim Sport LGT-232

Technical Specifications

	CE Version	FDA Version
Pulse Duration	50-300 μ s	200-400 μ s
Pulse Frequency	1-120Hz	1-120Hz
Intensity	0 - 80mA	0 - 100mA
Battery Life	More than 4 hours	More than 4 hours

Longest

Guangzhou Longest Science & Technology Co., L
Address: 5&6F, Building B4, No.11, Kaiyuan Avenue, Science City
Guangzhou Hi-tech Industrial Development Zone, 510530 Guan
Tel: 86-20-66353999 Fax: 86-20-66353920
E-mail: export@longest.cn Website: www.longest.cn



Health Life Longest Care

Longest



MStim Sport LGT-232

MStim Sport LGT-232 uses low frequency electric stimulation (NMES and TENS), help professional athletes or fitness, sports enthusiasts who with higher requirements on muscle performance, help them enhance muscle ability, improve athletic performance, relieve muscle soreness after exercise.

Applications

- ④ Body shaping
- ④ Relieve exercise pain
- ④ Increase athletic performance: endurance, strength, explosive force



Device Place

- ④ Electrodes Placement Guidance

A large number of muscle patch diagram to facilitate the correct placement of muscle position



Pocket for electrodes


Charger Socket

Four devices can be used and charged at the same time

Features

- ④ Portable device. Not only can do treatment with the app, but can also be used without app.
- ④ Wireless connection. Magnetic buckle electrode.
- ④ Charge Box

- ④ 7 Programs

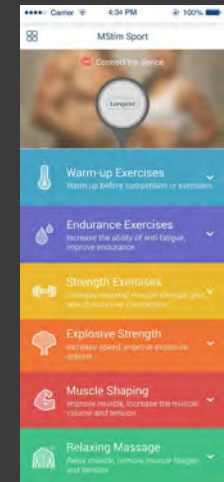

- ④ Intelligent muscle recognition

Different parameters to stimulate

Muscle contraction vibration

Sensor detection

Comfortable pulse width and intensity

With intelligent muscle recognition function, automatic output pulse width and intensity, no need to adjust any parameters, you can start muscle training.